

# Charging options for ReSound hearing aids

ReSound GN

ReSound offers rechargeable solutions with inductive charging and long battery life to keep your patients on the go.



Premium charger

Standard charger

Desktop charger

Custom charger\*

ReSound OMNIA™ Mini RIE	06	06	06	-
ReSound OMNIA™ RIE ReSound ONE™ RIE	01	01	01	-
ReSound OMNIA™ BTE ReSound ONE™ BTE	-	-	70	-
ReSound OMNIA™ Power/High Power BTE ReSound ONE™ Power/High Power BTE	-	-	80	-
ReSound OMNIA™ Customs Customs by ReSound	-	-	-	C
ReSound LiNX Quattro™ RIE ReSound Key™ RIE	03	-	03	-

\* Custom charger includes custom-made inserts. Remind patients to bring their custom-made charger along with their custom hearing aids to all follow-up appointments.

## How to match ReSound RIE and BTE hearing aids with the correct charger

Each charger type is offered with different inserts – known as the charging form factor (CFF) – which relate to a specific ReSound hearing aid collection or style.

The CFF is indicated with a number (see chart above) and is shown inside the charger and on the charging bay on the charger packaging.





# Charger features

	ReSound OMNIA			ReSound OMNIA ReSound ONE					ReSound OMNIA Customs by ReSound	ReSound LINX Quattro ReSound Key	
	Premium charger	Standard charger	Desktop charger	Premium charger	Standard charger	Desktop charger	Desktop charger	Desktop charger	Custom charger	Premium charger	Desktop charger
	Mini RIE	Mini RIE	Mini RIE	RIE	RIE	RIE	Standard BTE	Power/High Power BTE	Customs	RIE	RIE
Inductive charging	•	•	•	•	•	•	•	•	•	•	•
Portable with onboard power bank	•			•						•	
Portable without power bank		•			•				•		
USB-C connector and cable	•	•	•	•	•	•	•	•	•	•	•
LED on the back	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Lights up when the charger is connected to the main power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Shows when the charger is connected to a power supply
LED on the front	Shows the hearing aids charging			Shows the hearing aids charging					Shows the hearing aids charging	Shows the hearing aids charging	

## How to tell if the hearing aids are charging

Place the hearing aids in the charger and let the light indicator on the hearing aids guide you. For rechargeable Customs, the LED on the front of the charger will be your guidance.

You know that charging has started when the hearing aid lights blink slowly.

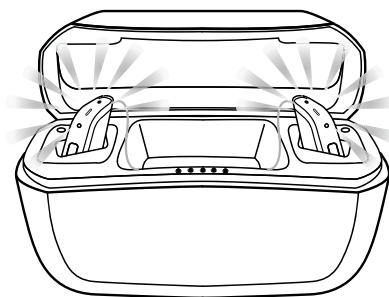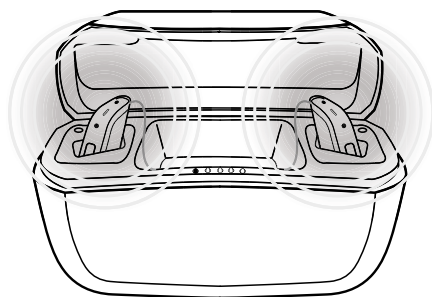

When fully charged, the hearing aid lights remain constant, not blinking.

The lights turn off when the hearing aids are removed from the charger.

Premium charger illustrated. The process is the same for the Standard and Desktop chargers. For more information, please see the individual user guide for each charger.

### Manufacturer according to FDA:

**GN ReSound North America**  
8001 E Bloomington Freeway  
Bloomington, MN 55420  
USA  
1-800-248-4327  
pro.resound.com